

Classic and Digital Slit Lamps

outstanding optics and bright, long-life LED lighting



bon

Slit Lamps SL-85V / SL-75V / SL-75 Elite

The right tool for every professional.



A slitlamp from bon - the right choice.

The excellent high resolution optics of our slit lamps have convinced thousands of professionals worldwide for decades.

The state of the art LED light source ensures perfect illumination for each and every examination and is energy saving. The LEDs promise a service life of up to 50,000 hours. They guarantee excellent colour rendition at any illumination level. The brightness control is comfortably integrated into the slit lamp base - full control at your fingertips.

The smooth mechanics of the slit lamp base allow you single handed control.

The diagnostic capability is enhanced by an integral yellow filter and the slit decoupling feature (SL-75V/Elite only).

DIN EN ISO 13485 certified quality

All models are produced in the EU by the world's biggest slit lamps manufacturer especially for bon. They have been engineered at the highest European standards.

Their ergonomic design, high quality optics and precision mechanical parts provide the user with an unparalleled experience while conducting a slit lamp examination.

These slit lamps can more than meet the demands of ophthalmic practice, even those tricky endothelium observations and without any additional light source.



C.S.O. S.r.l.,
Via degli Stagnacci 12/E,
50010, Badia a settimo (FI), Italy



bon SL-85 V

This slit lamp in Z-type construction (illumination from below) offers a bright LED illumination and a 5-fold magnification changer from 6x to 40x.

The illumination is controlled directly at the base of the unit.

The SL-85V is already prepared to take a DigiPro camera module. For this purpose, it offers a joystick trigger and R/L recognition in the base as standard.

bon SL-75 V

This model is the universal ophthalmic tool and is used by many ophthalmologists worldwide.

The SL-75V also offers modern LED illumination, a 5-fold magnification changer up to 40x and can be easily be upgraded to a digital slit lamp.

The classic H-type design with illumination from above allows the slit to be tilted horizontally by up to 20° (slit decoupling).



bon SL-75 Elite

The exclusive "Elite" slit lamps are the absolute top models of the bon slit lamp family. 40 years of experience in ophthalmology have resulted in the most modern premium class instruments.

This slit lamp convinces - in addition to the features of the SL-75V - with tangibly optimised mechanics with unique smoothness and smooth detents.

Its excellent, multi-coated optics guarantee the highest depth of field during the examination.



DigiPro makes quality visible.

Experience photos, videos and full-screen live images in unparalleled sharpness of detail. Even demanding endothelial and fundus shots are possible without the use of flash.
- Simply at the touch of a button!

The DigiPro has been the world's most successful digital slit lamp for many years, now in its fifth technology generation.


The SL-85V, SL-75V and SL-75 Elite slit lamps are already pre-equipped with a trigger in the joystick and automatic R/L detection as standard. Together with the camera module, they become a fully-fledged DigiPro in no time at all.

The camera has a professional 3/4 inch chip that has been specially developed for use in ophthalmology. It provides superior transmission values compared to a conventional beam splitter. An ambient illumination is included in the scope of delivery.

An automatic white balance ensures exceptionally true-colour image reproduction. The advanced „Selective Sharpness“ algorithm minimises imaging noise and improves details at high detail at high magnifications. The clever Multi-Shot function helps you to always capture the intended image.

DigiPro can be integrated into common practice applications via a professional DICOM interface.

A modern USB3.0 interface guarantees fast data transfer and high PC compatibility.

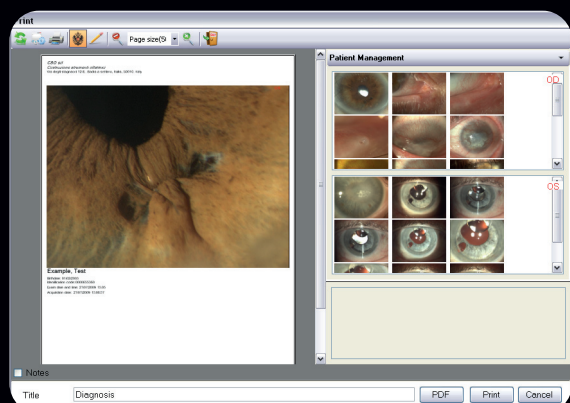

The Phoenix Software




The DigiPro camera module includes the advanced Phoenix software. Its network-compatible database is unlimited and ensures easy and professional image management.

The integrated image processing enables convenient marking, drawing, annotating and measuring directly in the picture. You can compare up to 6 images – optimal for case documentation.

Phoenix is a multi tasking software platform for the DigiPro slit lamp, the Cobra fundus camera, the Sirius Scheimpflug camera, the topographer Antares and the Perseus specular microscope – a world of diagnostic capability.

Technical Data SL-85V / SL-75V / SL-75 Elite

Microscope type	Galileian convergent with magnification change system
Eyepiece convergence angle	6°
Eyepieces	12,5x - suitable for spectacle wearers
Eyepiece adjustment	±8 D
Magnifications	6x 10x 16x 25x 40x
Corresponding real magnifications	5.6x 8.5x 14.8x 25.6x 39.3x
Field of view	41 mm to 5,7 mm
Interpupillary distance	50 mm to 80 mm
Integrated barrier filter	yellow
Voltage	15V DC 1A
Light source	White LED
Brightness	284000 LUX continuous adjustment
Color temperature	3100 K
Left/Right detection	integrated in instrument base
Trigger Button for digital camera	integrated in joystick
DigiPro camera module	PC connection USB3.0 Resolution 5MP (2448x2048 px) Frame rate 35 fps Pixel size 3.45 x 3.45 µm²

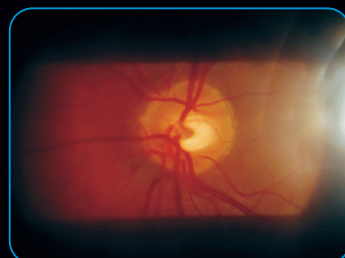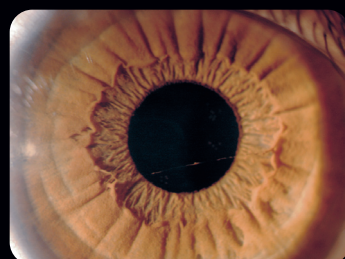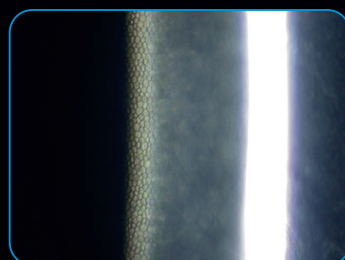

Slit Lamp SL-85V

Slit Lamps SL-75V and SL-75 Elite

Slit projection	1.3X	1X
Slit width (setting continuous)	0 – 15 mm	0 – 12 mm
Slit length (setting continuous)	1 – 15 mm	1 – 12 mm
Apertures	15, 9, 5.5, 0.3 mm	12, 9, 5, 3, 1, 0.2 mm
Filters	Blue, Red, Green (Red free)	Blue, Red, ND 50%, Green (Red free)
Slit rotation	± 90° continuous, on the TABO system	± 90° continuous, on the TABO system
Slit angle	0° horizontal	variable 0° 5° 10° 15° 20°
Rotation interval of the slit projection	±90°, angular scale, reference on 0°	±90°, angular scale, reference on 0° and ±10°
Working distance	69.5 mm	80 mm
Dimensions (HxWxD)	440 x 313 x 335 mm	675 x 313 x 335 mm
Weight (w/o camera)	7 kg	7.8 kg



Applanation Tonometer (optional accessory)

These precision measuring instruments are offering highest quality at a favorable bon price. You will feel no difference, neither in handling nor at accuracy. All bon tonometers are using standard measuring prisms, also available as disposables.

Polaris - The next Generation „TearScope“ (optional accessory)

As a practical addition to your slit lamp, Polaris extends your diagnostic possibilities in a simple way: With its LED cold light source, it enables you to perform interferometric analyses of the lipid layer of the tear film and to check the corneal surface. The Polaris can simply be attached to the slit lamp base.



Ergonomic Viewing Angle (optional accessory)

This eyepiece adapter is ideal for particularly tall users. It promises them relaxed working without neck pain. The viewing angle is 20°. The angle piece can be easily mounted on all current slit lamp models.

Rev. 2022-11-14

b o n Optic Vertriebsgesellschaft mbH

Stellmacherstr. 14 · 23556 Lübeck, Germany · Phone 49 451 80 9000
Fax +49 451 80 900-10 · E-Mail call@bon.de · Internet www.bon.de

bon